

Connect Clip

TIPS & TRICKS



Oticon's Connect Clip is a multifunctional accessory, which can serve as a microphone, intermediary streaming device, or remote control. Consider each feature to make the most of your Connect Clip.

PAIRING

For each of these functions, the Connect Clip must first be paired to the hearing aids.

1. Turn the hearing aids off and on. The hearing aids will be in "pairing mode" for 3 minutes.
2. Turn the Connect Clip on by holding the multi-button for 3 seconds. The light will blink **yellow**, while searching for the hearing aids.
3. Place the hearing aids close to the Connect Clip (no more than 8 inches apart) and wait.
4. A steady **yellow** light on the Connect Clip confirms a successful pairing.



REMOTE MICROPHONE

Use the Connect Clip to send the speaker's voice directly to the hearing aids.

Consider this function when listening from a distance (up to 65 feet) or in noise (in the car, at the park, during walks).

1. After pairing the hearing aids, press and hold the multi-button for 3 seconds. The Connect Clip light will turn solid **magenta** to confirm "Remote Mic Mode."
2. Speaker should wear the Connect Clip approximately 8 inches below his or her mouth. Ensure the microphone is not covered by clothes or jewelry.
3. To mute/unmute the microphone, hold the multi-button for 3 seconds.




STREAM AUDIO FROM ANDROID


Use the Connect Clip to stream audio (music, video, phone call) to the hearing aids.

Clip the Connect Clip to your shirt and use the microphone for hands-free phone calls.

1. After pairing the Connect Clip to the hearing aids, the Connect Clip must be paired to the Android device.
2. Press the multi-button and the volume-up button. Hold for 6 seconds until the Connect Clip light blinks **blue**.
3. Enable Bluetooth on the Android device and search for "ConnectClip".
4. Select "ConnectClip". If the Android requests a code, enter "0000". The light will turn solid **blue** to confirm "Bluetooth Mode."
5. Once paired to an Android phone, calls will automatically stream to the hearing aids, and the indicator light will change to **green**. Use the Connect Clip microphone for hands-free calls.
6. Music and other audio will automatically stream to the hearing aids, and the indicator light will change to **orange**.

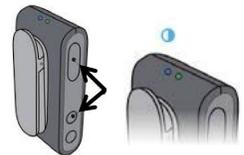

REMOTE CONTROL

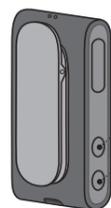

- **Volume Up**
- short press
- **Volume Down**
- short press



- **Next Program**
- Press and hold for 3 secs for next program
- **Mute hearing aid microphones - long press**
short press Volume Up or Down will unmute



www.cdhe.isdh.in.gov